

INSTRUCTION FOR USE

SENSIS*Strip* Pecan nut 20/5 Tests (Cat. nr. HU0030129/HU0030169)



Lateral-flow Device for the Determination of Pecan nut in Food and as Cleaning Control Monitoring

Sensitivity for food matrix	2.5 ppm
Sensitivity for swabbing	0.007 µg/cm ²
Sensitivity for rinse water	0.17 mg/L

1. GENERAL INFORMATION

The pecan tree (*Carya illinoensis*) belongs to the family of Juglandiaceae. The amount of protein in pecan nut is about 9%. Some of these proteins are known for being allergenic, like the 2S albumin Car I 1 or the 11S legumin Car I 4. Compared to other nuts the allergenic potential of pecan nuts is still slightly characterized. Pecan nut allergies are relatively seldom, but can be very distinct in particular cases. For pecan nut allergic persons hidden pecan nut allergens in food are a critical problem. Already very low amounts of pecan nut can cause allergic reactions, which may lead to anaphylactic shock in severe cases. Because of this, pecan nut-allergic persons should strictly avoid the consumption of pecan nut containing food. Cross-contamination, mostly in consequence of the production process, is often noticed. This explains why in many cases the existence of pecan nut residues in food cannot be excluded. For this reason, sensitive detection systems for pecan nut residues in food-stuff are required.

The **SENSIS*Strip* Pecan Nut Lateral Flow Device** represents a sensitive detection system and is particularly capable to detect pecan nut residues in food matrices, rinse water and swabs.

2. PRINCIPLE OF THE TEST

The **SENSIS*Strip* Pecan nut** test is based on the principle of immunoassay. Pecan nut containing sample is given into a reactions vial containing biotinylated antibody directed against pecan nut proteins. After 3 minutes incubation at room temperature a test strip is placed into the reaction vial. The sample migrates along the nitrocellulose membrane by capillary forces. Along its way it releases gold nanoparticles conjugated to streptavidin. An antibody-gold complex is

formed. For positive samples a red line is formed when the liquid reaches the test line area. In case of negative samples, no line is formed. In any case, above the test line area a red control line appears, indicating the validity of the test. The test is evaluated after another 5 minutes.

3. PRECAUTIONS

Full compliance of the following good laboratory practices (GLP) will determine the reliability of the results:

- 1) Store the kit at 2-8°C.
- 2) Do not use the kit after its expiry date.
- 3) Prior to beginning the assay procedure, bring all samples and reagents to room temperature (20-25°C).
- 4) Extraction buffer should be mixed by gentle inversion or swirling prior to use. Do not induce foaming.
- 5) Once the assay has been started, all subsequent steps should be completed without interruption and within the recommended time limits.
- 6) Replace caps in all the reagents and samples immediately after use.
- 7) Use separate disposable consumables for each transfer of sample to the reaction vial in order to prevent cross-contamination.
- 8) Do not mix components from different batches.
- 9) Do not use reagents after expiration date.

NOTE: The swab sampling device included in this kit may be supplied as sterile with a sterility expiration date printed on the device. However, this kit does not require a sterile sampling device, there-fore the swab sterility expiration date does not affect the kit expiration date and can be disregarded.

4. CONTENTS OF THE KIT

The kit contains components and reagents for 20 tests or 5 tests. They have to be stored at 2-8°C. Expiry data are printed on the labels of the reagent containers and the outer package.

Content	20-strip	5-strip
Test Strips, in tube with desiccant stopper	20 pcs	5 pcs
Reaction vials	20 pcs	5 pcs

Extraction tubes with caps	20 pcs	5 pcs
Extraction Buffer, 60 mL, ready-to-use.	1 pcs	1 pcs
Disposable Pipettes, 0.3 mL	20 pcs	5 pcs
Disposable Pipette, 3 mL	1 pcs	1 pcs
Disposable Spatulas	20 pcs	5 pcs
Swab Sticks	20 pcs	5 pcs
Evaluation Card	1 pcs	1 pcs
Tubes and vials racks	by kit box	by kit box
IfU	1 pcs	1 pcs
QR-Code for evaluation with RapidScan ST5 lateral flow strip reader	1 pcs	1 pcs

5. SAMPLE PREPARATION

Due to high risk of cross-contamination all applied instruments like applicator, mortar, vials etc. have to be **cleaned thoroughly** before and after each sample. Allergen proteins adhere very strongly to different surfaces. In certain cases, they can resist a common dishwasher cleaning. To identify possible cross-contamination caused by previous extractions it is strongly recommended to note the sequence of the extractions for pattern recognition.

Chocolate and other products with high polyphenol content tend to show reduced results. To overcome this effect a special extraction additive can be ordered separately (HU0030100).

5.1 Solid samples / Liquid samples

- 1) Homogenize sample using appropriate methods depending on its specific nature (e.g. grind, crush, mix).
- 2) *Solid samples:* Transfer one spatula of sample to an extraction tube. Alternatively, in order to increase precision, weigh out 0.2 g of sample into an extraction tube.

Liquid samples: Transfer a half spatula of sample liquid to the extraction tube. Alternatively, in order to increase precision, pipette 0.2 mL of sample into an extraction.

- 3) Add 3 mL of ready-to-use extraction buffer to the sample by using the disposable 3 mL pipette.
- 4) Close extraction tube with cap and shake for 1 minute.
- 5) Let the solid remains sediment. Depending on nature of the samples this might take 1-2 minutes. Alternatively centrifuge at 2000 g or higher.
- 6) Remove cap and transfer 0.3 mL of sample supernatant into a reaction vial using a disposable 0.3 mL pipette.


5.2 Rinse water

- 1) Adjust the pH of the sample to 7 (+/- 0.5).
- 2) Transfer 0.3 mL of sample into a reaction vial using a disposable 0.3 mL pipette.

5.3 Swabbing samples

DRY SURFACES

- 1) Mark out 5x5 cm area or use swab directly on (e.g. uneven) area.
- 2) Transfer 1 mL of ready-to-use extraction solution into an extraction tube by using the disposable 3 mL pipette.
- 3) Moisten a swab by dipping into the tube.
- 4) Swab marked area by using crosshatch (1. horizontally, 2. vertically, 3. diagonally) technique while rotating the tip.

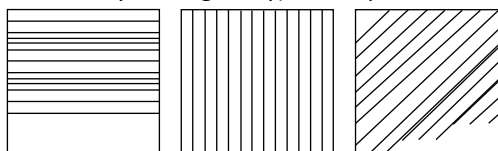

- 5) Place swab into the tube and break off the tip.
- 6) Close extraction tube with cap and shake for 1 minute to release the sample from the swab.
- 7) Remove cap and transfer 0.3 mL of sample supernatant into a reaction vial using a disposable 0.3 mL pipette.

WET SURFACES


Apply same method as described for dry surfaces without prior need to moisten the swab.

6. PROCEDURE

- 1) Prepare samples as described above.
- 2) After transfer of the sample to the incubation vial add cap and shake for 15 seconds. Make sure that the biotinylated antibody is completely dissolved.
- 3) Incubate for 3 minutes.
- 4) Remove cap and place one strip into the vial. For proper strip orientation make sure that the arrows on the cover foil point downwards.
- 5) Incubate for 5 minutes.
- 6) Remove strip from the vial and evaluate immediately.

7. EVALUATION

SENSIS*Strip* lateral-flow devices are evaluated according to the following scheme:

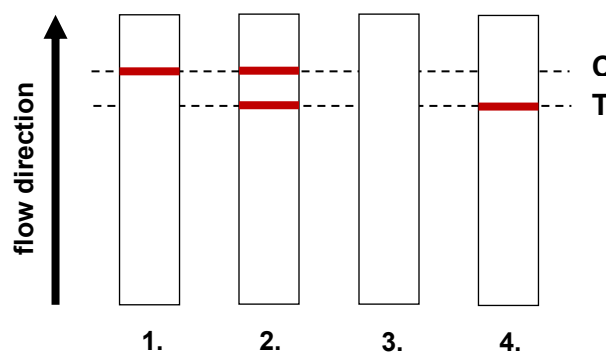

1. **Negative:** visible control (C) line, no test (T) line

2. **Positive:** visible control (C) and test (T) lines
3. **Invalid:** neither control (C) and test (T) lines visible
4. **Invalid:** no control (C) line and visible test (T) line

For a better distinguishing between negative, borderline and positive samples a colour card for evaluation is provided with the kit. The intensity of the test line has to be compared with the different increments of the colour card. Results lower than increment 3 should be treated as negative. Results according increment 3 or higher should be treated as positive. Since the increments of the colour card are ranging up to 10 a semi-quantitative evaluation is also possible. This can be improved by taking into account the results stated in the validation report of the product.

8. PERFORMANCE

8.1 Sensitivity

LOD (total Pecan nut) of the SENSIS*Strip* lateral-flow test is 2.5 ppm for food matrix, 0.17 mg/L for rinse water and 0.007 µg/cm² for swab samples applying the procedure above. The corresponding amounts of Pecan nut protein can be calculated by anticipating a protein content of Pecan nut of 9%.

NOTE: Sensitivity may vary depending on matrix and processing of a complex food mixture. For achieving reliable results each matrix should be validated prior to routine testing.

8.2 Cross-reactivity

For the following foods no cross-reactivity could be detected:

Adzuki bean	Curcuma	Oyster
Almond	Dill	Paprika
Apricot	Duck	Pea
Barley	Egg white	Peach
Bean, white	Egg yolk	Peanut
Bovine	Ewe's milk	Pepper
Bovine gelatine	Fennel	Pine nut
Buckwheat	Fenugreek	Pistachio
Caraway	Flaxseed	Poppy seed
Cardamom	Garden cress	Pork
Carob bean	Garlic	Potato
Carrot	Gliadin	Pumpkin seed
Cashew	Goat's milk	Radish
Cayenne	Guar gum	Rice
Celery	Hazelnut	Rye
Cherry	Horseradish	Sesame
Chestnut	Kidney bean	Shrimp
Chia	Kiwi	Soy flour
Chicken	Lamb	Soy lecithin
Chickpea	Leek	Soy milk

Chili	Lentil	Split peas
Cinnamon	Lupin	Sucrose
Clove	Macadamia	Sunflower seed
Cocoa	Milk powder	Thyme
Coconut	Mustard, yellow	Tomato
Cod	Nutmeg	Turkey
Corn	Oats	Wheat
Cow's milk	Onion	White cabbage
Cumin		

The following cross-reactions were determined:

Brazil nut	0.125%
Walnut	100%

8.3 High-dose-hook Effect

Reduced or absent signals can occur in case of very high concentrations. The test gives valid results up to a concentration of 10000 ppm for food samples, according 26.7 µg/cm² for swabs and 667 mg/L for rinse water samples.

8.4 Additional Performance Data

Additional data can be found in the corresponding validation report of the product, which can be inquired at Gold Standard Diagnostics Budapest.

9. LIABILITY

Gold Standard Diagnostics Budapest shall not be liable for any damages to the customer caused by the improper use of the kit and for any action undertaken as a consequence of results.

Gold Standard Diagnostics Budapest shall not be liable for the unsafe use of the kit out of the current European safety regulations.

Distributed by Emport LLC
 Interested in auto-shipments?
 Other products? Contact us!
orders@emportllc.com
 412-447-1888 or 866-509-4482